

GenePasQ Data Sheet / Reading guide

Title _____

DATE ()

Name ()

• Gene amplification by PCR


OneStep RT-PCR kit

(Lot #) _____


①	50 °C	30 min	} _____ 39 cycles
②	95 °C	15 min	
③	95 °C	15 sec	
④	60 °C	45 sec	
⑤	4 °C		

• GenePasQ

(Lot #) _____

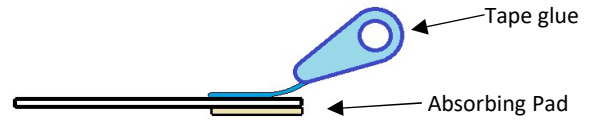

Paste the GenePasQ strips which is glued behind of absorbing pad with Tape glue.

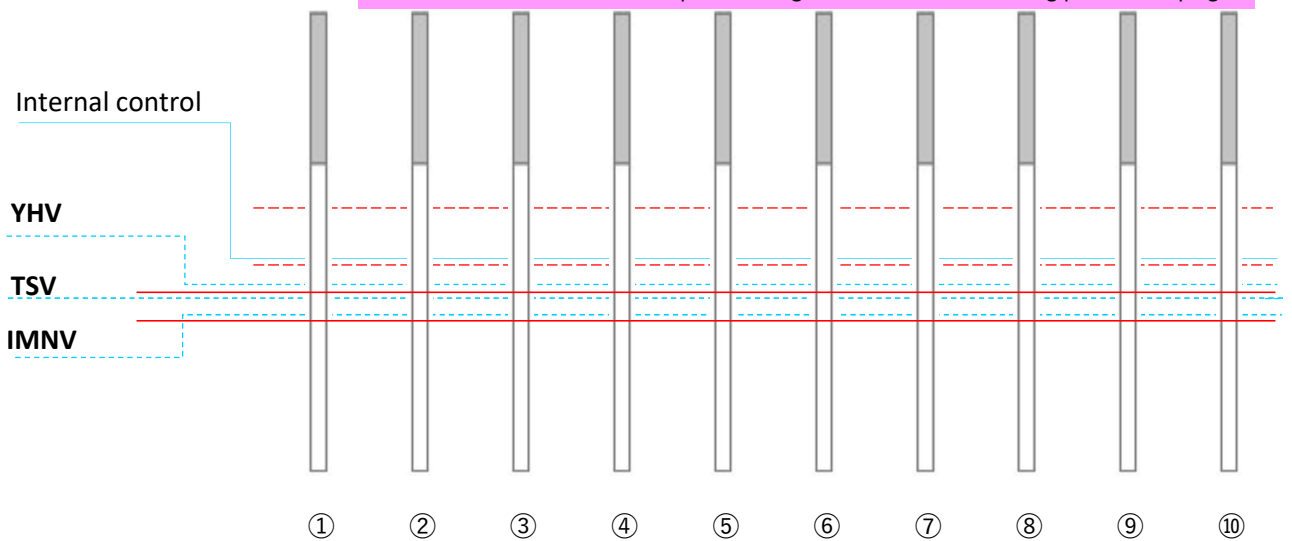

	①	②	③	④	⑤	⑥	⑦	⑧	⑨	⑩
YHV										
TSV										
IMNV										

* Internal control (IC) line often becomes faint when a target gene is amplified.

test	comment
①	
②	
③	
④	
⑤	
⑥	
⑦	
⑧	
⑨	
⑩	